

explorē 8

User Guide



Table of Contents

Chapter 1 – Overview.....	4
Introduction	4
Box Contents	5
Additional Kit Contents	5
Physical Description	6
Using the Folding Stand (Kit option).....	9
Getting Started.....	11
Chapter 2 – Using Your System.....	13
Powering On	13
Powering Off	13
Magnifier Mode	14
Zooming In and Out	14
Changing Contrast.....	14
Turning lights on and off	15
Lines and Blinds	15

Autofocus	16
Freezing an Image.....	16
Navigating on a Frozen Image.....	17
Gallery Mode	18
Accessing the Gallery	18
Deleting Images from Gallery	19
Using an External Screen	19
Chapter 3 – Exporting Images to Computer	20
Chapter 4 – Menu.....	21
Accessing the Menu and Changing Settings	21
Chapter 5 – Updating Your System.....	26
Chapter 6 – Recharging Your System.....	26
Chapter 7 – Troubleshooting.....	29
Appendix A – Safety and Maintenance	31
Appendix B – Specifications.....	34

Appendix C – FCC Information 35

Appendix D – Warranty 37

Chapter 1 – Overview

Introduction

Thank you for purchasing the HumanWare explorē 8, the light and easy to use portable electronic handheld magnifier with a large, easy to read 8” screen. With your explorē 8, you’ll be reading those fine prints in a matter of seconds.

Being magnifier users ourselves, we at HumanWare know that in addition to the best available image quality, a good magnifier needs to be simple, portable, comfortable to use, and should always be ready when you need it, where you need it. Your explorē 8 was designed with just that in mind, for we take great pride in offering you a quality product that will meet both your needs and expectations. Carefully read the operating and safety instructions described in this user guide prior to using your device. If you have any questions about your system, please contact HumanWare Customer Service by consulting

the contact information in the back of this user guide.

Copyright 2024 by Technologies HumanWare Inc., all rights reserved. No part of this publication may be reproduced or distributed in any form or by any means, or stored in a database or retrieval system, without the prior written permission of Technologies HumanWare Inc.

Box Contents


- explorē 8 electronic magnifier
- User Guide
- USB charger
- USB cable
- High-Resolution video cable
- Protective case
- LCD protector
- Microfiber cloth

Additional Kit Contents


- explorē 8 folding stand
- carry bag

Physical Description


Front face

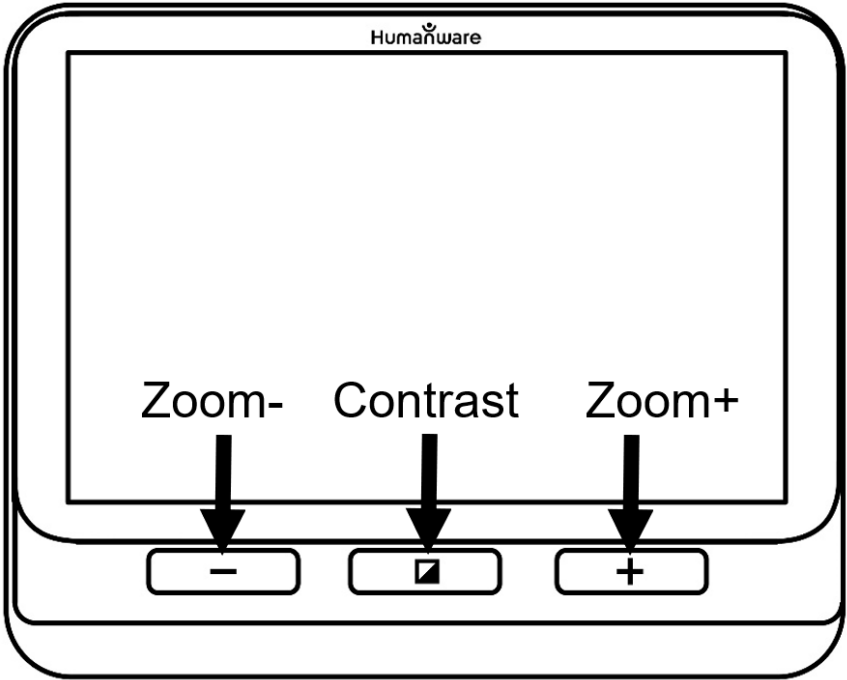


Left side

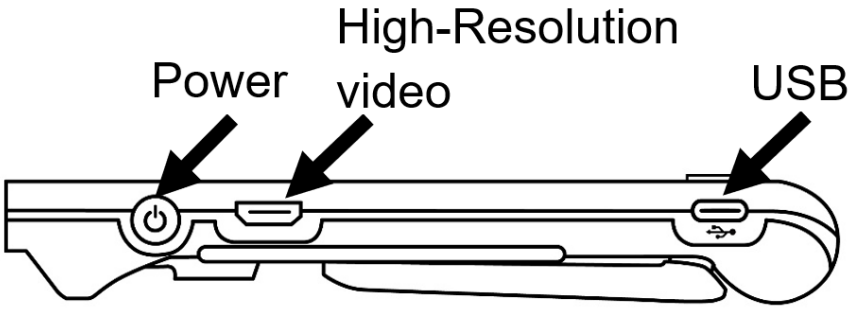

Top edge

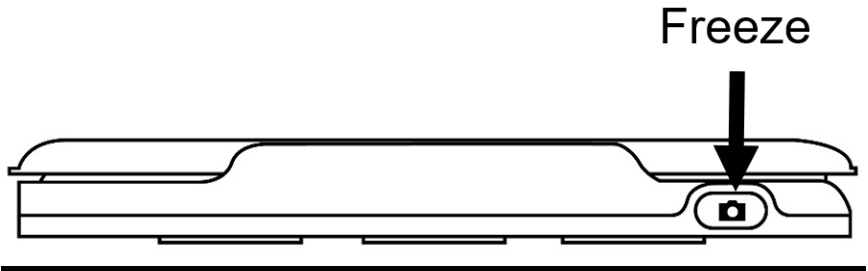

Back side

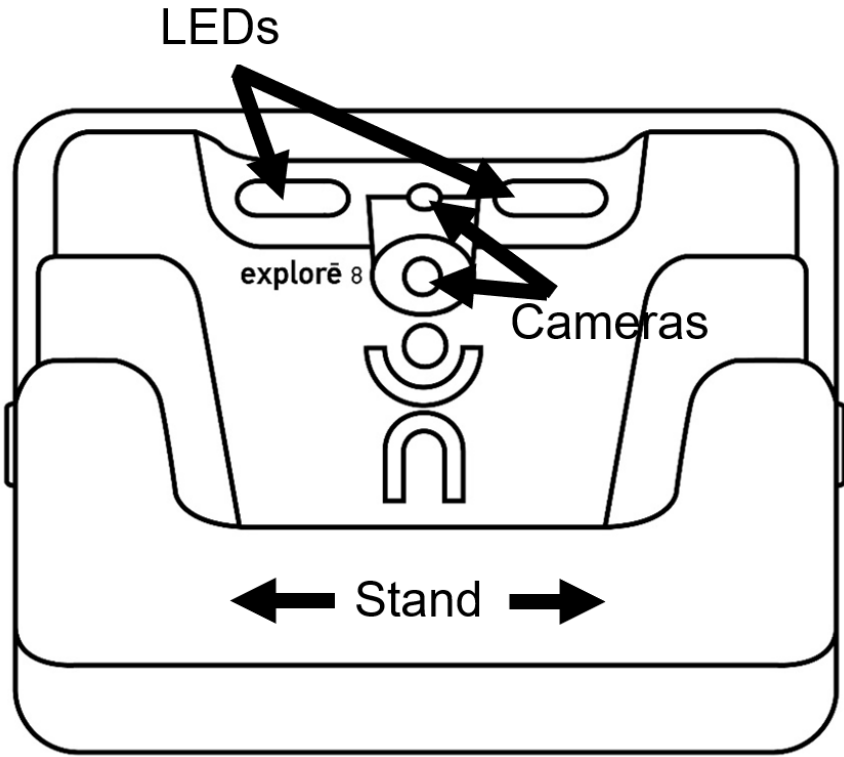

Using the Folding Stand (Kit option)


When used with its foldable stand, **explorē 8** rises above the table, allowing for plenty of writing room.

Unfold the (optional) folding stand, then slip the device into the stand cradle to use it in "Stand Mode." (see Folding stand images 1 and 2)


When you are finished working in "Stand Mode," slide the device out from the cradle using the reverse motion.

To optimize "Stand Mode," refer to the settings adjustments found in [Section 4](#).

1.

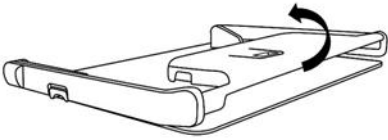


2.

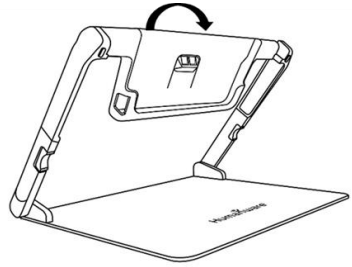

3.




4.

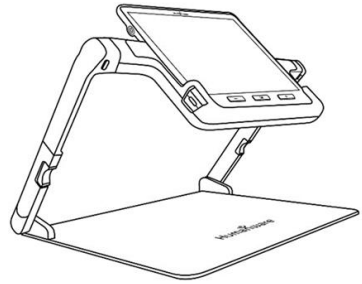


5.




6.

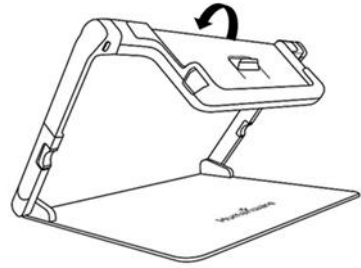

7.

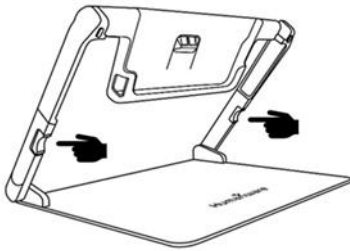

8.

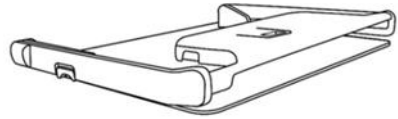

Getting Started

1. Connect the USB cable to the USB charger.
2. Hold explorē 8 so that the screen side is facing you. Carefully connect the cable to

the USB port located on the left side of your device and plug the charger into an outlet.

IMPORTANT: Use caution while inserting the cable into the USB port.


3. Let your device recharge until your battery is full. The device will become warmer during the charge – this is normal.
4. Once the device is fully charged, carefully disconnect the cable.
5. Press **Power** or unfold the reading stand to start using the device.

Chapter 2 – Using Your System

Powering On

To start your magnifier, simply press **Power** on the left side of the device. The device can then be held by the sides, just like you would hold a tablet.

Alternatively, you can turn on your device by simply unfolding the reading stand behind the device and placing it over what you wish to read.

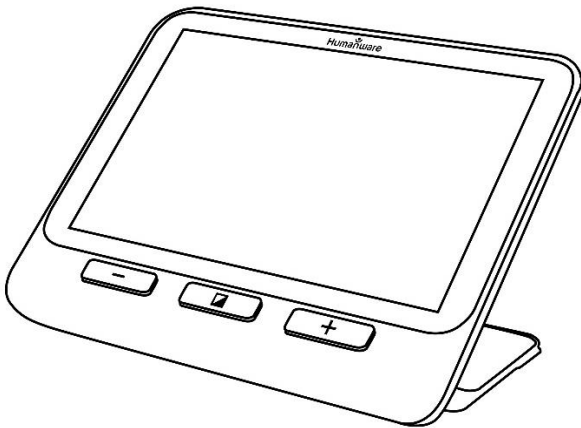

After a few seconds, the explorē 8 screen will turn on and display the image from the camera. This is the magnifier mode.

Powering Off

To turn off your device, simply press **Power**.

Magnifier Mode

Magnifier mode is explorē 8's main mode of operation. With the stand closed, you will be able to magnify distant objects. With the stand open, you will be able to magnify close objects or documents on a table.

Zooming In and Out

To zoom in, simply press or press and hold **Zoom+** until the image is of the appropriate size. Similarly, you can zoom out by pressing or pressing and holding **Zoom-** until you are satisfied with the image size.

Alternatively, you can enable the **pinch and zoom** gesture (decrease and increase the distance between two fingers) from the menu to zoom in and out of a live image (see Ch. 4).

Changing Contrast

explorē 8 allows you to customize contrast to better suit your needs and preferences. Simply press **Contrast** to change contrast until you find one you are comfortable with.

The available contrast modes are listed below:

- Full color, Black on white, White on black, Black on yellow, Yellow on black

Other contrast modes are also available from the menu (see Ch. 4):

- Black on orange, Orange on black, Black on cyan, Cyan on black, Black on purple, Purple on black, Black on green, Green on black, Blue on yellow, Yellow on blue, White on blue, Blue on white, Red on white, White on red, Red on black, Black on red

Turning lights on and off

Press and hold **Zoom+** and **Zoom-** simultaneously to turn off the lights if necessary, and to turn them back on.

Lines and Blinds

With the explorē 8, you can add a horizontal line or blinds over your image for easier reading.

To add or remove a line, press and hold **Contrast** and **Zoom+** simultaneously.

To add or remove blinds, press and hold **Contrast** and **Zoom-** simultaneously.

Autofocus

explorē 8 features an autofocus that can be set to **locked**, **unlocked**, and **manual** from the menu (see Ch. 4).

When autofocus is set to **unlocked** (by default), explorē 8 will dynamically adjust its focus to the object currently being viewed. The focus can be forced by touching anywhere in the screen.

When autofocus is set to **locked**, or you are using your magnifier in tabletop position, explorē 8 will maintain its current focus regardless of any movements or changes.

When autofocus is set to **manual**, you can touch any position on the screen to have explorē 8 focus on that particular location. Note that the **manual** setting only works in distance viewing mode. Tabletop mode will continue to be locked when **manual** mode is selected.

Freezing an Image

explorē 8 enables you to freeze an image while in magnifier mode for easier reading, an especially useful feature when reading labels

on products on stores' top shelves. To freeze an image, simply press **Freeze** while aligning your magnifier towards what you wish to read. This will take a snapshot of what your device was oriented towards. The snapshot will be displayed on your screen until you decide to return to Magnifier mode.

To return to magnifier mode, simply press **Freeze** once again.

Navigating on a Frozen Image

To zoom in on still images, simply press or press and hold **Zoom+** until the image is of the appropriate size. Similarly, you can zoom out by pressing or pressing and holding **Zoom-** until you are satisfied with the image size.

Alternatively, you can enable the **pinch and zoom** gesture (decrease and increase the distance between two fingers) from the menu to zoom in and out of a frozen image (see Ch. 4).

Hold your finger down and slide it across the screen to see other parts of a frozen image.

As for a live image, press **Contrast** to change contrast until you find one you are comfortable with.

Gallery Mode

To save a live or frozen image to your Gallery, simply press and hold **Freeze** for a half-second while in Magnifier mode. You should see a Save message displayed briefly on your screen meaning your image is now saved inside explorē 8's Gallery.

Alternatively, you can enable the **tap and hold** gesture (keep your finger down on the screen) from the menu to save an image in the Gallery (see Ch. 4).

Accessing the Gallery

To access your images in the Gallery, press **Contrast** and **Freeze** simultaneously. Swipe your finger left or right on the screen to select an image. You can also press the arrows on the screen to select an image.

You can zoom in and out of saved images by pressing or pressing and holding **Zoom+** and **Zoom-**. Alternatively, you can enable the **pinch and zoom** gesture (decrease and increase the distance between two fingers) from the menu to zoom in and out of a saved image (see Ch. 4).

To exit the Gallery and return to Magnifier mode, press **Back** (on screen) or press and hold **Contrast**.

Deleting Images from Gallery

To delete an image from the Gallery, access the Gallery and navigate to the image you wish to delete (by swiping your finger left or right on the screen or using the arrow buttons), then press and hold **Freeze**.

Alternatively, you can enable the **tap and hold** gesture (keep your finger down on the screen) from the menu to delete an image (see Ch. 4).

You will have the option to Delete or Cancel. Press **Delete** to confirm deletion. To cancel deletion, press **Cancel** or **Contrast**.

You can delete all your images at once from the menu (see Ch. 4). Press and hold **Contrast** and select **Delete Gallery**.

Using an External Screen

For a larger image, you may want explorē 8's content to be displayed on an external screen (like a TV). To do so, simply connect your explorē 8 to an external screen using the included High-Resolution video cable. Please

refer to Appendix B for a list of supported High-Resolution video formats.

Chapter 3 – Exporting Images to Computer

To export Gallery images saved on explorē 8 to your computer, connect your device to your computer with the USB cable provided. You should see a “PC connected” message appear on your device. You can now copy and paste the images from your Gallery to your computer using your computer’s file browser.

When connected to a personal computer, the device will act exactly like a USB thumb drive or an external drive. Gallery images are located in the Images folder on your explorē 8.

When you are copying or deleting images from the device, never disconnect the USB cable as some images could be lost in the process.

Chapter 4 – Menu

Accessing the Menu and Changing Settings

Press and hold **Contrast** to enter the menu. You can navigate the menu by swiping up or down with your finger on the screen, pressing the up/down arrows on the screen, or pressing **Zoom+** to scroll up and **Zoom-** to scroll down. You can also press the up and down arrows on the screen.

To select an item, press on it. To go back to the previous menu, press **Contrast** or **Back**. Press **Contrast** again to return to Magnifier mode. To quickly exit the menu, press and hold **Contrast** or **Back**.

The menu contains the following items:

- **Gallery** access Gallery images
- **Autofocus** Unlocked*, Locked (Locked* in stand mode), Manual (Locked in stand mode)
- **Screen Brightness** (Change the screen's brightness level): 1 to 9 (8*)
- **Touchscreen** (submenu, tap to enable or disable gestures and vibrations)

- **Pinch & zoom** (On*)
(decrease/increase distance between two fingers to zoom out/in)
- **Vibrations** (On*)
- **Tap & hold** (On*) (hold finger down in Magnifier to save image to Gallery)
- **Temporary focus lock (toggle)** (On*)
Available only when using the optional folding stand. When turned on, double-tap on the screen with one finger to lock focus in its current state. This feature prevents the device from refocusing, allowing tasks like document signing without continuous refocusing of the camera. Double-tap a second time with one finger to deactivate.
- **Image overview (toggle)** (Off*)
Usable only with the optional folding stand. The two-finger “Overview” gesture is available in document viewing mode: Place two fingers near the top of the screen and slide them to the bottom while maintaining screen

contact. The view will zoom out to 1x. This allows you to orient the magnified area on the screen. Remove both fingers from the screen to return to the previous zoom level.

- **Toggle button (toggle)** (Off*) Must be set to On to use the explorē 8 on its optional folding stand. Use the “T” (text) icon on the folding stand for viewing items while using the stand or tap the mountain icon and remove the device for handheld or distance viewing (icons are displayed at the top left of the screen).
- **Auto-off** (Choose how long the device stays on): 1 minute, 2 minutes, 5 minutes*, 10 minutes, Never
- **Preferred Zooms** (Users can select the default zoom level when the device is powered up):
 - **Base:** Previous, 2x* to 30x
 - **Folding Stand:** Previous, 2x* to 15x
 - **Distance:** Previous, 2x* to 30x

- **Contrast**

- **Preferred Contrast** (select contrast in which your explorē 8 starts up): Previous* (last used contrast), Color, enabled contrasts
- **Menu Colors** (change menu colors): Black on white*, White on black, Black on yellow, Yellow on black, Black on orange, Orange on black, Black on cyan, Cyan on black, Black on purple, Purple on black, Black on green, Green on black, Blue on yellow, Yellow on blue, White on blue, Blue on white, Red on white, White on red, Red on black, Black on red
- **Enabled Contrasts** (select enabled color pairs for magnifier mode): Black on white*, White on black, Black on yellow, Yellow on black, Black on orange, Orange on black, Black on cyan, Cyan on black, Black on purple, Purple on black, Black on green, Green on black, Blue on yellow, Yellow on blue, White on blue, Blue on white, Red on white, White on red, Red on black, Black on red

- **Lines and blinds:** Line, Blinds, None*
- **Lights:**
 - **Base:** On*, Low, Off
 - **Folding Stand:** On*, Low, Off
 - **Handheld:** On, Low, Off*
- **Menu Text Size:** Small, Medium*, Large
- **Language:** English US*, English UK, French (CA), French (FR), Spanish (US), Spanish (ES), German, Dutch, Danish, Polish, Portuguese (PT), Portuguese (BR), Norwegian, Italian, Swedish, Finnish, Russian, Czech, Japanese, Korean, Icelandic, Estonian, Turkish, Lithuanian, Hungarian, Latvian, Arabic, Chinese (Simplified), Chinese (Traditional), Hebrew, Georgian, Farsi, Kurdish
- **Display:** 50Hz , 60Hz*
- **Delete Gallery**
- **Factory Defaults**
- **About** (View information about your device): Date and time, Software, Hardware, Serial Number

Chapter 5 – Updating Your System

First, go to:

<https://www.humanware.com/exploresupport>

Then, download the latest version of the explorē 8 software. Connect your device to your computer using the provided USB cable. Using your computer's file explorer, paste your downloaded file directly into your explorē 8 root. When the transfer is completed, unplug your USB cable from your computer. Restart your explorē 8. Your device should automatically update itself.

Chapter 6 – Recharging Your System

IMPORTANT: The battery becomes warm during recharging. This is normal. Do not recharge the battery near a heat source.

To recharge the battery using a household power outlet:

- Connect the USB cable to explorē 8's USB charger. Then, carefully connect the

USB cable to your device and plug the charger to an outlet.

- Let your device charge until your battery is full. It should take approximately 3.5 hours to recharge an empty device.

During the recharge process, the image below will show momentarily when the battery starts charging or when the device is turned on.



The lightning bolt icon means your device is being charged. The number to the right represents the percentage of battery charge remaining.

Due to lower voltage, charging from the computer USB may be slower. In some cases, your device may not charge at all depending on the amount of power supplied by the computer. The most effective way to recharge the device is to use the supplied USB cable

connected to the USB charger and recharge it from a household power outlet.

Although it may lengthen the charging time, the device can be used when it is recharging.

Please note that if the battery is fully discharged or hasn't been used for a very long period, there will be a delay of up to a few minutes before the device shows any activity. This is normal.

Also note that like all electronic devices, your explorē 8 will become warm when in use or when being charged – this is normal.

Battery safety precautions:

- Do not disassemble or modify the battery.
- Use only the specified HumanWare charger.

There is a risk of overheating, fire or explosion if the battery is put in a fire, heated, subjected to impact, put in contact with water, or if its terminals are shorted.

Chapter 7 – Troubleshooting

explorē 8 does not power on:

Press and hold **Power**. If the unit still does not turn on, your device's battery may need recharging.

Connect the provided USB cable with the USB charger into a power outlet or just the USB cable into a running computer. If the battery was completely discharged, it may take a few minutes before the device shows any signs of activity (see Ch. 6 for more info).

Screen is black:

If the unit is laying on a table or any flat surface, the image may be black. Lift your magnifier to see whether this is the reason why your screen is black.

If this does not solve the problem, try restarting your device.

Image is out of focus in magnifier:

Make sure Autofocus is set to Unlocked in the menu. Move the device back and forth over the document you are viewing. This will make explorē 8 re-focus on the document.

explorē 8 shuts down after inactivity:

This is a normal behaviour intended to preserve the battery. By default, explorē 8 will shut down after 5 minutes of inactivity. It is possible to configure the default auto-off timer to 1 minute, 2 minutes, 5 minutes or 10 minutes.

Colors are difficult to detect:

First make sure there is sufficient ambient light around you.

Choose a different color combination by pressing **Contrast** while in magnifier mode.

There is a lot of glare:

Tilt the object you are looking at in a different direction to reduce glare. You can also try

turning off the lights by pressing and holding **Zoom+** and **Zoom-** simultaneously.

explorē 8 does not respond:

Make sure explorē 8 is not connected to power via USB cable; if it is, unplug it. Press **Power** for 10 seconds. The system will shut down. Then, press Power for one second to turn explorē 8 back on.

Other problems:

Restore to factory defaults in the menu.

Press and hold **Contrast**. Scroll down to Factory Defaults and select Yes.

Appendix A – Safety and Maintenance

You must not stare at the LED light source located on the back side of the magnifier.

Keep away from water or high humidity. Do not submerge.

Power off your explorē 8 before cleaning it.

Use only cables provided with your explorē 8.

Battery:

The explorē 8 has an internal rechargeable Li-Ion battery. The battery must be replaced by a qualified HumanWare technician.

- Do not use or leave the unit near heat source and fire. Do not store in high temperatures.
- Only use USB-compliant power supplies to charge the battery.
- Do not disassemble or modify the unit.
- Do not immerse the unit in water, do not allow it to get wet.
- Do not strike or throw the unit.
- Do not pierce or strike the unit with sharp objects or a hammer.
- If the battery leaks and liquid gets into the eyes, do not rub eyes. Instead, rinse the eyes with clean running water and seek immediate medical attention to prevent injury.
- If the battery leaks and comes into contact with skin, immediately rinse area with clean running water to prevent injury.

Cleaning the LCD screen:

Make sure your device is off and unplugged.

Make sure you only use microfiber cloths to clean the LCD screen like the one included with your device.

Carefully and gently wipe the microfiber cloth across the LCD screen.

If you choose to clean your LCD screen with a liquid cleaner, make sure this cleaner is appropriate for LCD screens. Distilled water, isopropyl alcohol and a vinegar-water solutions (no more than 10% vinegar) are safe cleaning liquids for your LCD screen.

Use only the minimal amount of liquid to prevent spills inside the device. If liquid seeps in, damages to the device might occur. Make sure no liquids remain on your screen. Repeat the previous steps for spots that are harder to remove.

Never use household cleaners to clean your explorē 8.

Never spray any liquid directly on your LCD screen.

Never use cleaning products that contain ethyl alcohol, ethyl acid, ammonia, acetone or methyl chloride.

Never use products that contain isopropyl alcohol to clean the lenses of the cameras.

Never use paper towel to clean your LCD screen.

Appendix B – Specifications

Size: 24 x 204 x 153.2 mm

Display: 8" ISP LCD

Weight: 580 g

Power Input: USB Type C connector.
5V@3A, 9V@2A or 12V@1.5A

Battery life: 3.5 hours

Charge Time: 2 hours (when device is off)

High-Resolution video output

Compatibility:

Supported High-Resolution video output formats:

1280x720p 60Hz, 50Hz

1920x1080p 60Hz, 50Hz, 24Hz

1920x1080i 60Hz, 50Hz

No support of DVI emulation

NOTE: *Hi-resolution performance may vary from specifications with the use of some external screens or monitors.*

Operating conditions:

Temperature: 10 °C to 40 °C / 50 °F to 104 °F

Humidity: 20% to 80%

Transport and storage conditions:

Temperature: -20 °C to 60 °C / -4 °F to 140 °F

Humidity: 10% to 95%

Appendix C – FCC Information

FCC Statement:

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Warning:

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.

—Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

—Consult the dealer or an experienced radio/TV technician for help.

Note: This unit was tested with shielded cables on the peripheral devices. Shielded cables must be used with the unit to ensure compliance.

Appendix D – Warranty

HumanWare warrants explorē 8, effective from the date of purchase, to be free of any defects in material and workmanship for a period of 2 years unless otherwise required by law in the country or region of purchase.

This warranty is non-transferable and applies to all cases where the damage is not a result of improper use, mistreatment, or negligence. Improper use is use of this device other than described in this manual. In no event shall HumanWare or its distributors be liable for indirect or consequential damages.

No replacement or repair covered by the warranty will be carried out unless the system is accompanied by a copy of the original bill of purchase. Please keep your original receipt. If the system has to be returned, please use the original packaging.

This page intentionally blank.

This page intentionally blank.

This page intentionally blank.

This page intentionally blank.

Customer support:

For customer support, please contact the HumanWare office nearest you or view our Website at: www.humanware.com

North America / Amérique du Nord:

1 (800) 722-3393

support@humanware.com

Europe:

(0044) 1933 415 800

eu.support@humanware.com

Australia / Asia - Australie / Asie:

(02) 9686 2600

au.sales@humanware.com